



# ETG 255



## APPLICATION INSTRUCTIONS



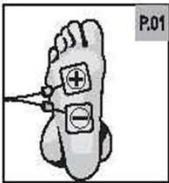
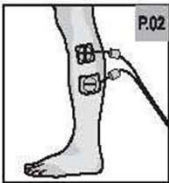
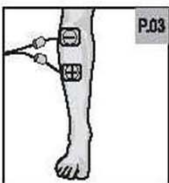
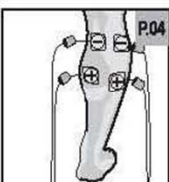

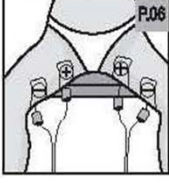
- EMS/TENS Device ETG 255
- Electric muscle stimulation through EMS application
- Pain alleviation through stimulation current/TENS
- 2 channels with adjustable intensity
- LED display for each channel
- 15 user programs for muscle stimulation
- 3 user programs for nerve stimulation
- Beginner, advanced and professional mode
- Timer function and display background illumination



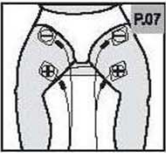
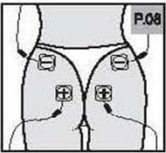
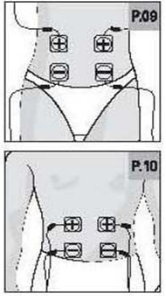
## APPLICATION EXAMPLES FOR ADHESIVE ELECTRODES

### Stimulation positions

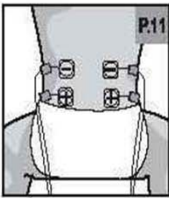
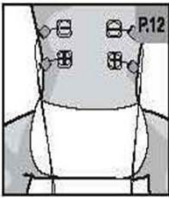
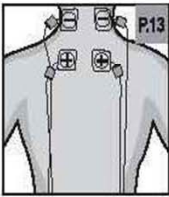
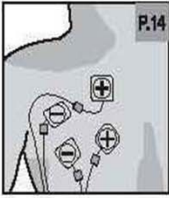
The stimulation position depends on the muscle groups to be stimulated. For the different suggested positions, refer to the pictograms next to the images for electrode positioning. The following table lists the different muscle groups as well as some useful information about the best stimulation position and how to intentionally cause a contraction (muscle tension).

| Muscle group                     | Placement of electrodes   | Stimulation positions   | Manually causing a contraction (muscle tension)  |
|----------------------------------|---|---|--|
| Muscles on the soles of the feet |    | Sitting position, place feet on the floor   | Forcefully tension the muscles in the soles of your feet and try to dig your feet into the floor   |
| Fibula muscles                   |    | Sitting position, place feet on the floor   | Powerfully tension the fibula muscles by pressing your big toe onto the floor and lift the outer toes off the floor at the same time           |
| Front Tibia muscles              |  | Sitting position, place your feet beneath a piece of furniture so that the ankles can no longer be bent | Powerfully tighten the front fibula muscles by pressing the tips of the feet against an obstacle that counteracts this movement. Press upwards |
| Calf muscles                     |  | Sitting position, so that back and feet are supported. You best sit in a doorframe                      | Powerfully tighten the calf muscles by pressing the tips of the feet against an obstacle that counteracts this movement.                       |
| Rear thigh muscles               |  | Lie flat on your stomach, ankles fixed into place without feeling uncomfortable                         | Powerfully tension the rear thigh muscles by trying to bend your knees   |
| Muscles for pulling up your legs |  | Sitting position, place a hard object between the knees without feeling uncomfortable                   | Powerfully tension the adductor muscles by trying to press your knees together   |

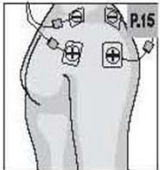
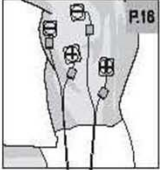
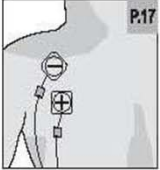
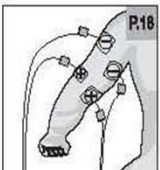
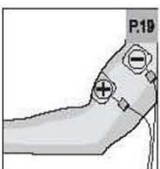
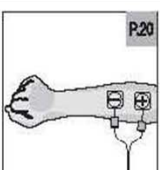
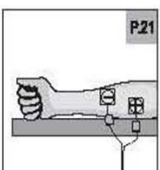
## APPLICATION EXAMPLES FOR ADHESIVE ELECTRODES

| Muscle group        | Placement of electrodes  | Stimulation positions  | Manually causing a contraction (muscle tension)  |
|---------------------|--|--|--|
| Front thigh muscles |   | Sitting position. There are two variants for this exercise: static, block knee movement, or dynamic, move gains an obstacle, use heavy weights for this purpose  | Powerfully tension the front thigh muscles by trying to stretch your legs                                    |
| Gluteal muscles     |   | Lie on your stomach or stand up.   | Powerfully tighten your gluteal muscles by contracting them and trying to push your thighs behind your torso |
| Stomach muscles     |  | Lie on your back (you may lift it slightly). There are two variants for this exercise: static, for this purpose simply stat muscle contraction by performing the motion described on the side; or dynamic by additionally moving the torso towards the thighs; in this case, make sure not to place the focus on the lumbar spine (backward curvature); the knees should always be firmly pressed together | Tension your stomach muscles by trying to lift your head and shoulders off the floor                         |

## APPLICATION EXAMPLES FOR ADHESIVE ELECTRODES

| Muscle group               | Placement of electrodes   | Stimulation positions   | Manually causing a contraction (muscle tension)  |
|----------------------------|---|---|--|
| Lower back muscles         |    | <p>Sitting position.<br/>Please note: Due to the anatomic peculiarity of the lower back muscles, training in this mode requires very strong muscles.<br/>Place the electrodes in the region of the back muscles as illustrated.</p> | <p>Powerfully tighten the lower back muscles by trying to sit as upright as possible</p> |
| Back muscles               |    | <p>Sitting position.</p>  | <p>Powerfully tighten the lower back muscles by trying to sit as upright as possible</p> |
| Cervical vertebrae muscles |   | <p>Sitzende Position.</p>   | <p>Powerfully tighten the lower back muscles by trying to sit as upright as possible</p> |
| Trapeze muscles            |  | <p>Sitting position.</p>  | <p>Tighten the trapeze muscle by trying to powerfully raise and lower your shoulders</p> |

## APPLICATION EXAMPLES FOR ADHESIVE ELECTRODES

| Muscle group            | Placement of electrodes   | Stimulation positions   | Manually causing a contraction ( muscle tension)  |
|-------------------------|---|---|---|
| Shoulder joint muscles  |    | Sitting position, elbows in the arm rest so that the arm rest poses an obstacle when moving away from the body  | Powerfully tighten the shoulder muscles by pressing the elbows away from the body.                      |
| Large back muscle       |    | Sitting position, elbows outside the arm rest so that the arm rest poses an obstacle when moving towards the body   | Powerfully tighten the large back muscle by pressing the elbows towards the body.                       |
| Chest muscles           |    | Sitting position.<br>The palms of the hand touch.<br>Warning notice concerning the placement of electrodes in the thorax (coronary) region: increased risk of triggering cardiac fibrillation | Powerfully tighten your chest muscles by pressing the palms of your hands together                      |
| Rear Upper arm muscles  |   | Sitting position, place hands and forearms on a table.  | Powerfully tighten the rear upper arm muscles by pressing the palms of your hands on the table.         |
| Front Upper arm muscles |  | Sitting position, forearms are placed on a table, the palms must face upwards. Fix your elbows so that they do not move during the stimulation.   | Powerfully tighten the front upper arm muscles by moving the palms of your hands towards your shoulders |
| Hand stretchers         |  | Sitting position, forearms and palms are placed on the table.   | Powerfully tighten the hand muscles by trying to lift your hands.                                       |
| Hand benders            |  | Sitting position, forearms are placed on the table<br>Take a durable, hard object into your hand while keeping your fingers slightly bent.  | Powerfully tighten the hand muscles by gripping the object in your hand                                 |